

TRUEresult Test Strips

available on prescription

diabetes
SIMPLE BLOOD GLUCOSE MONITORING

7CM
ANYTHING ELSE IS MORE OF A HANDFUL



a *Small* revolution
twist-and-go convenience for people with diabetes



© 11/12 NIPRO DIAGNOSTICS LTD TR11105B 11/12

<p>TRUEresult twist True On-the-Go Convenience</p> <p>PIP Code: Blue Energy - 353-3940 Pink Aura - 353-3957 Phantom Black - 353-3932</p>	<p>TRUEresult True Performance State-of-the-Art Design</p> <p>PIP Code - 353-3924</p>	<p>TRUEresult test strips State-of-the-Art Technology Delivers Precision and Accuracy</p> <p>Available on prescription PIP Code - 353-3973</p>
---	--	---

For accurate and precise results recommend the no coding, TRUEresult and TRUEresult twist Systems today!

Nipro Diagnostics (UK) Limited
25 Barnes Wallis Road, Segensworth East
Fareham, Hampshire, PO15 5TT
free upgrade ☎ 0800 08 588 08



in partnership



in partnership



TRUEresult
twist



a *Small* revolution
for accurate convenient testing for people on-the-go

Easy to Use – Just Twist on and Test Fast

- Just twist the **TRUE**result twist on to top of a new vial of strips
- Test with results in as fast as 4 seconds



Convenience Features Provide Competitive Advantages


- **No Coding** – eliminates a step that can lead to inaccurate results
- **Results As Fast As 4 Seconds** – means hardly any waiting time
- **Tiny, 0.5 Microlitre Sample Size** – less blood means less painful testing
- **Easy-to-Read Display** – for easy viewing of test results
- **Strip Release Button** – eliminates handling of strips after testing
- **Test Memory** – stores up to 99 tests



TRUEresult Test Strips

state-of-the-art technology ensures accurate results and promotes first-test success

Easy to Use – Just Twist on and Test Fast

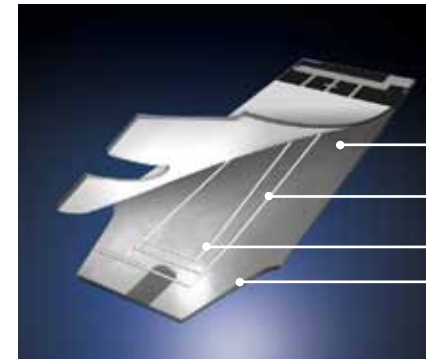

- Quad Electrode Laser Accuracy provides enhanced precision and accuracy of test results
- Two fill-detect electrodes ensure adequate blood sample, preventing false starts and erroneous results
- Two electrodes accurately measure the glucose reaction and send signal to sensor
- Innovative TRUEfill beveled tip allows for greater sampling precision

TRUEresult Test Strips are utilized by both TRUEresult and TRUEresult twist Systems, and feature the latest in testing technology

